




# CARE SURGICAL

## CS Universal Half-Roll

Reusable half-roll foam support designed for safe and comfortable positioning of the patient to reduce risk of pressure-related injury e.g. to support above the hollow of the patient's knee when in supine position.

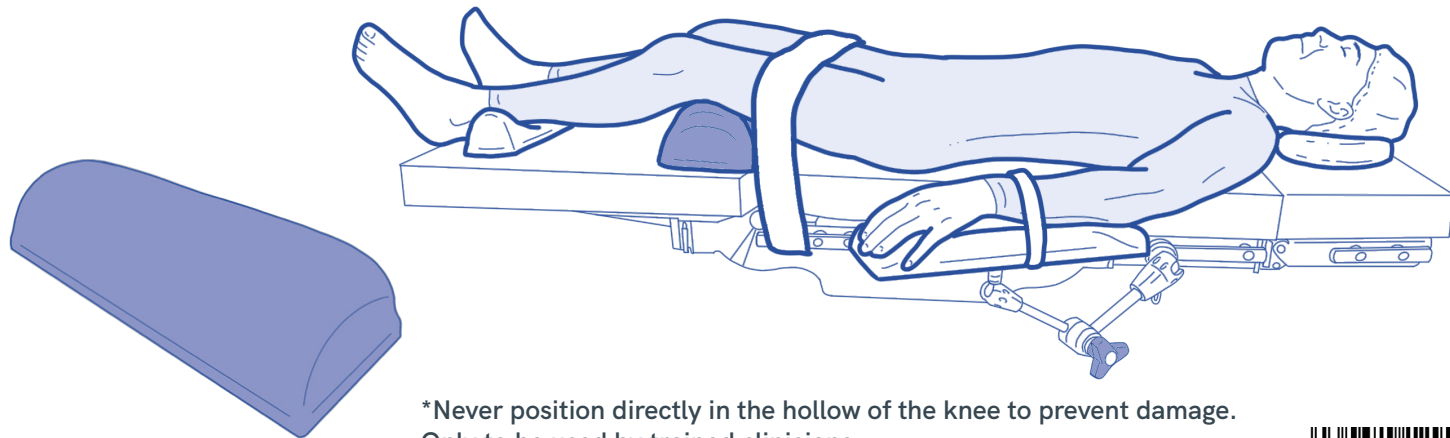

\*Never position directly in the hollow of the knee to prevent damage.  
Only to be used by trained clinicians.



Manufacturer

UK: Care Surgical Ltd  
Unit 6 Ringtail Road  
Burscough, L40 8JY  
Tel: 01704 336671

USA: Care Surgical LLC  
43 Broad Street  
Suite B106, Hudson,  
MA 01749, USA  
Tel: 866-243-4107

[www.care-surgical.com](http://www.care-surgical.com)

Care Surgical Supplies Limited, Office 2, 12a Lower Main Street,  
Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



UK Conformity  
Assessed



Reuse



CE Mark



Latex Free



Non Sterile



LOT Number



EC REP

Product Code: CSM-2907 Rev.4 02/24



## Head-to-Toe Check list

This check should be repeated continuously – at the latest when the position of the table has been changed. At all stages, check that no tubing is trapped between the patient and the supports and that all bony protrusions are padded.

### The legs and feet


- + Avoid compression of the back of the knee
- + Position a half roll above the back of the knee instead
- + No pressure on the calves to avoid compartment syndrome
- + Heels, fibular and tibial nerves should be completely free of pressure as well

### Arm adjacent to the body

- + The medial elbow joint should not lie on an edge
- + It should be placed against the thigh in a neutral position
- + Check again for intraoperative body movements

### The head

- + Head to be positioned centrally
- + Without rotation of the cervical spine
- + Head to be well padded (minimum pressure on the back of the head)

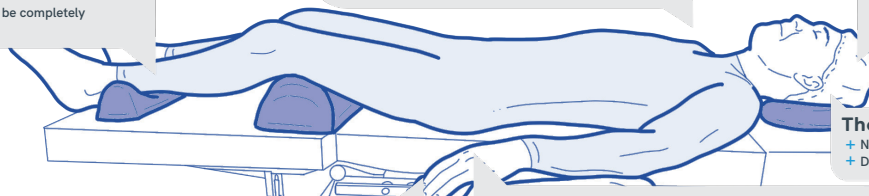

### The arms

- + Positioning of the arms should be done at the end
- + If possible, the arms should be out of position

### Arm abducted from the body

- + Fix the arm supports at shoulder level to the table
- + Avoid drop hand
- + Position upper arms above table level
- + With an abduction in the shoulder joint of 0-40°, the hand should be positioned in pronation, with 40-90° in supination.
- + An abduction of more than 90° should be avoided at all costs to prevent possible damage to the brachial plexus.
- + Arms should always be secured

### The neck:

- + Neutral position
- + Do not overextend

To see our full range of Care Surgical Medical Products, visit [www.care-surgical.com](http://www.care-surgical.com)



Manufacturer

UK: Care Surgical Ltd  
Unit 6 Ringtail Road  
Burscough, L40 8JY  
Tel: 01704 336671

USA: Care Surgical LLC  
43 Broad Street  
Suite B106, Hudson,  
MA 01749, USA  
Tel: 866-243-4107



Medical Device



UK Conformity  
Assessed



Reuse



CE Mark



Latex Free



Non Sterile



5 060624 129070

[www.care-surgical.com](http://www.care-surgical.com)



Care Surgical Supplies Limited, Office 2, 12a Lower Main Street,  
Lucan Co. Dublin, K78 X5P8 Ireland