




CARE SURGICAL

Product Code: CSM2930

CS Carbon Spine Extension- Merivaara™


⚠ Only to be used by trained clinicians. Ensure product is secure before positioning patient.

Fit the CS Carbon Spine Extension to the Merivaara Grand Promerix Operating Table to achieve imaging and patient positioning capabilities comparable to that of a specialist dual-rail spine table.

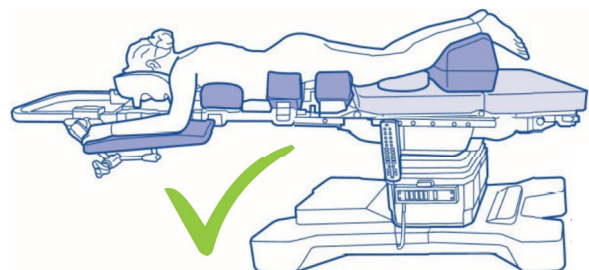


CSM-2930 CS Carbon Spine Extension- Merivaara™
Contains: CSM-2742, CSM-2730 & CSM-2737

CSM-2730 CS Carbon Spine Extension Frame

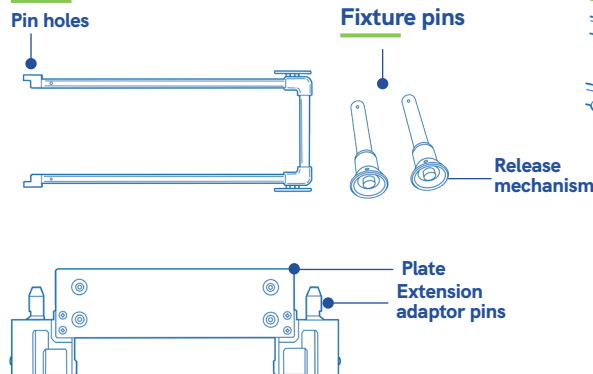

CSM-2737 CS Carbon Spine Extension Leg

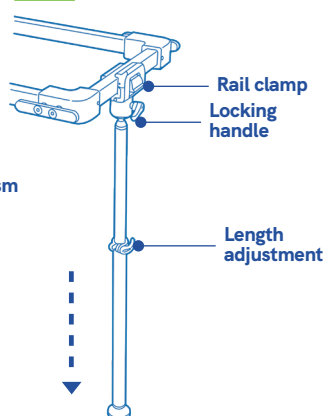


CSM-2742 CS Carbon Spine Extension Adaptor - Merivaara




Only to be used by trained clinicians. Recommended to be attached to the HEAD END of the operating table only for FULL WEIGHT BEARING CAPACITY. Only to be used with the operating table in the NEUTRAL position. Do not exceed the safe working load of 250kg on the HEAD END.


WARNING: If attachment at the FOOT END is required the SWL is REDUCED to 180kg and MUST be used with the stabilizing leg.

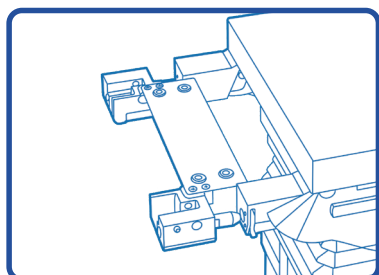


1 Line up the pins of the extension adapter with the operating table.

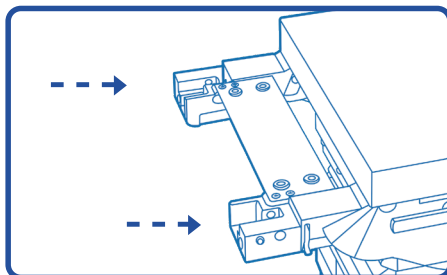


2 Insert the extension adapter into the operating table until the adapter locks into place. Ensure the extension adapter is flush against the operating table.

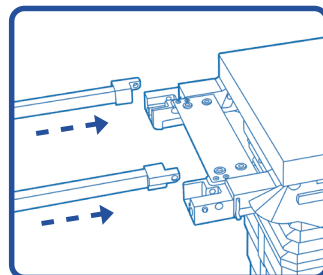

3 Supporting the full weight of the frame, line the pin holes of the extension rails up with the holes of the extension adapter.

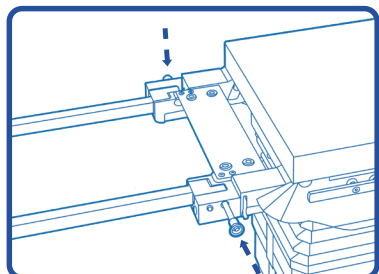

4 Insert the fixture pin- taking care to press the release mechanism in the centre of the pin in- through the extension adapter and extension rail to secure the extension.

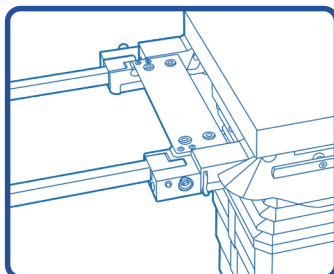

5 Repeat on opposite side and ensure both sides are secured. Release the weight of the frame into position.

Cleaning Instructions:



Clean using standard hospital wipes



Do not immerse in water



Do not autoclave



UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
Tel: 01704 336671

USA: Care Surgical LLC
43 Broad Street
Suite B106, Hudson,
MA 01749, USA
Tel: 866-243-4107



Maximum Patient Weight

www.care-surgical.com

Care Surgical Supplies Limited Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland

⚠ **SAFE DISPOSAL:** customers should adhere to any local laws and regulations to the safe disposal of medical devices.
⚠ **DO NOT USE** products that appear to be damaged or worn. In an event of product damage or patient injury, please contact Care Surgical.






CARE SURGICAL

Product Code: CSM2930


CS Carbon Spine Extension-Merivaara™

To adjust the leg:


Intended use: The CS Carbon Spine Extension Leg is intended to be used to support and stabilise the CS Carbon Spine Extension Frame during the procedure.

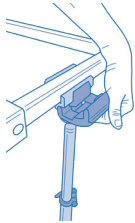


- 1** Place the leg in the centre of the head end of the CS Carbon Spine Extension Frame.

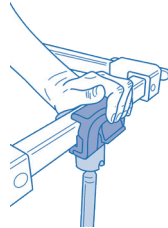

- 2** Ensure that the clip is closed and the leg is firmly attached to the frame.

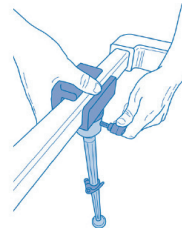





- 3** Ensure the locking handle is locked and secured by tightening clockwise.

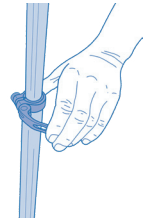

- 4** If the operating table needs to be raised, raise using the table hand control and then loosen the blue lever on the rod to set the new height before tightening the blue lever again.

WARNING: If the table needs to be lowered, first raise the table to remove any pressure on the leg. Then make sure that the black lever (fig. 4) is released before lowering the table. Once in the final position, tighten the black lever again.

WARNING: Do NOT remove the leg whilst the patient is in position without first raising the table to release any tension or weight on the leg. When the support leg has any load bearing, the lever (fig. 4) must not be released to prevent the patient from making a sudden downward movement. The only time the leg doesn't need to be engaged is when you want to maneuver the operating table.

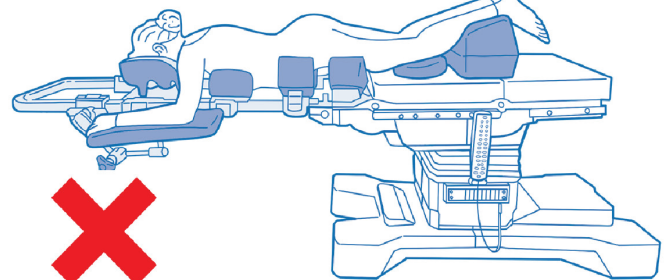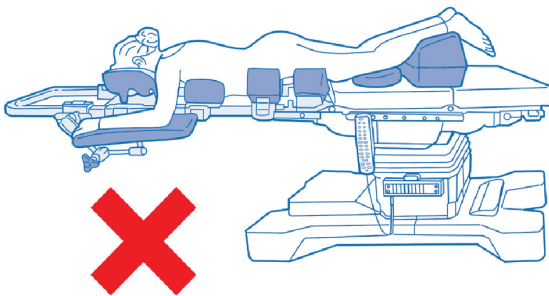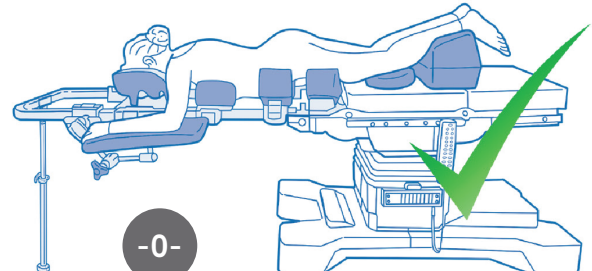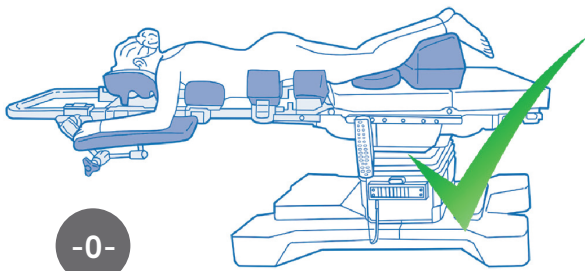

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
Tel: 01704 336671

USA: Care Surgical LLC
43 Broad Street
Suite B106, Hudson,
MA 01749, USA
Tel: 866-243-4107

www.care-surgical.com

Care Surgical Supplies Limited Office 2, 12a Lower Main Street,
Lucan Co. Dublin, K78 X5P8 Ireland



SAFE DISPOSAL: customers should adhere to any local laws and regulations to the safe disposal of medical devices.
DO NOT USE products that appear to be damaged or worn.
In an event of product damage or patient injury, please contact Care Surgical.

