




# CARE SURGICAL

# CS Egg Box Foam 80" x 20"

Reduces the chance of pressure related injury to the patient. Tear along the perforations for smaller pieces to fit the patient's arms, legs and chest.

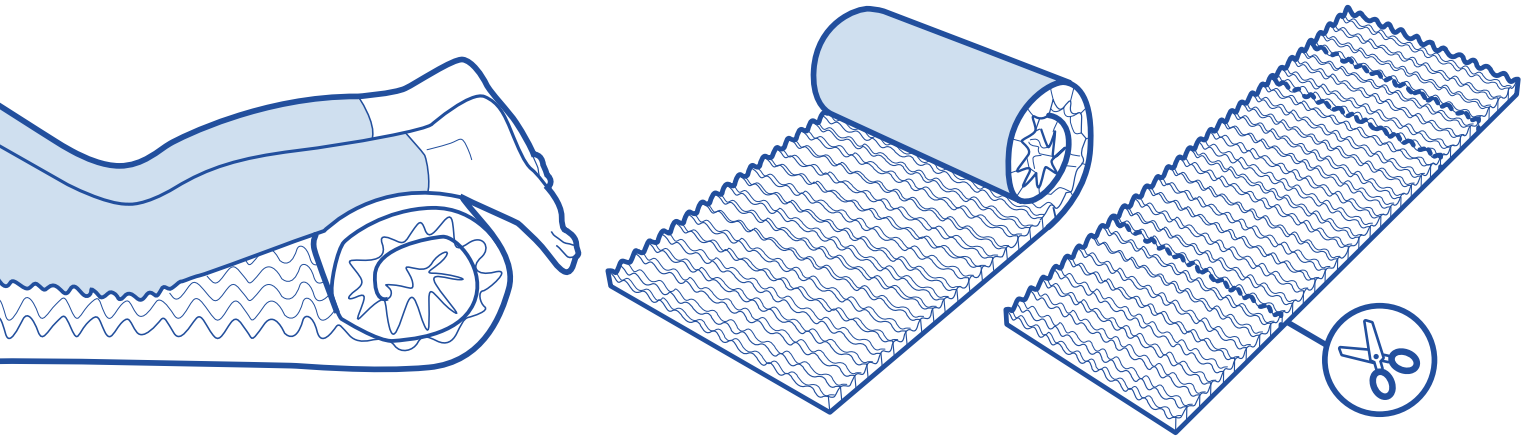

Sheet dimensions 80" x 20" / 200cm x 50cm



Manufacturer

UK: Care Surgical Ltd  
Unit 6 Ringtail Road  
Burscough, L40 8JY  
Tel: 01704 336671

USA: Care Surgical LLC  
Marlborough  
MA 01752  
Tel: 866-243-4107

[www.care-surgical.com](http://www.care-surgical.com)

Care Surgical Supplies Limited, 69 Esker Woods Drive  
Lucan, Co Dublin, K78PX45, Ireland



Do not reuse



CE Mark



Latex Free



Non Sterile



5 060624 120190

Australian Sponsor, Summit Surgical  
PO Box 60, Landsborough  
QLD 4550 Australia

LOT Number:

Product Code: CSM-2585 Rev 05 08/21

## GTIN

- + 5060624120190

## Description

- + Corrugated polyurethane foam sheet 80" x 20" x 2"

## Intended use

- + Padding for patients

## Indications for Use

- + Single Use
- + Can be used on operating tables, trolleys beds or chairs
- + Each sheet is 200x50cm
- + Can be cut to size, or tear along perforations
- + Reduces and distributes pressure to reduce the risk of pressure ulcers
- + The undulating structure of the foam enables air to circulate to assist in maintaining the patient's temperature
- + Latex and phthalate free
- + Instructions for use are packaged with the product

## Warnings, Precautions, Contraindications & Adverse Effects

- + Incorrect placement of this product could cause pressure marks or ulcers on the patient

## Packaging

- + Each unit is packaged in a recyclable polyethene bag
- + 1 unit is packaged into a cardboard box for sale

## Transport & Storage

- + Transport and store in dry conditions and out of direct light where possible
- + Exposure to light may cause yellowing of areas of the product. This is purely cosmetic and will not affect the quality or the safety of the product