




CARE SURGICAL

CS Lateral Head Support

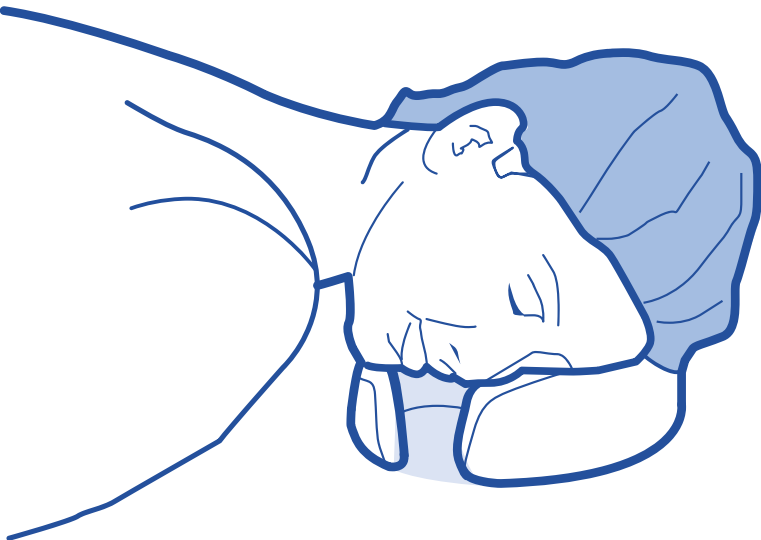

It is recommended that the patient's head is turned every 2 to 4 hours per ICS guidelines

The CS Lateral Head Support is a soft, wipeable lateral head support for patient's in the lateral position.

No need to disconnect the tubing.



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
Tel: 01704 336671

USA: Care Surgical LLC
Marlborough
MA 01752
Tel: 866-243-4107

www.care-surgical.com

Care Surgical Supplies Limited, 69 Esker Woods Drive
Lucan, Co Dublin, K78PX45, Ireland

EC REP



Do not reuse



CE Mark



Latex Free



Non Sterile



5 060624 120404

Australian Sponsor, Summit Surgical
PO Box 60, Landsborough
QLD 4550 Australia

LOT Number:

Product Code: CSM-2812 Rev 02 06/21



CARE SURGICAL

CS Lateral Head Support



Cleaning instructions:



Clean using standard hospital wipes



Do not immerse in water



Do not autoclave



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
Tel: 01704 336671

USA: Care Surgical LLC
Marlborough
MA 01752
Tel: 866-243-4107

www.care-surgical.com

Care Surgical Supplies Limited, 69 Esker Woods Drive
Lucan, Co Dublin, K78PX45, Ireland



Do not reuse



CE Mark



Latex Free



Non Sterile

Australian Sponsor, Summit Surgical
PO Box 60, Landsborough
QLD 4550 Australia



5 060624 120404

Product Code: CSM-2812 Rev 02 06/21

GTIN

- + 5060624120404

Description

- + Horseshoe-shaped head support made from polyurethane foam

Intended use

- + Supports patient's heads in the lateral position

Indications for Use

- + Reusable - clean using non-alcohol hospital wipes
- + Recommended for adults, but use will depend on patient's stature
- + Supports patient's head in the lateral position
- + Can be used for patients who are in the prone or lateral position
- + Reduces the incidence of pressure ulcers and marking
- + Patient can be repositioned without disconnecting their tubing, so PEEP pressure is maintained
- + Latex and phthalate free
- + Instructions for use are packaged with the product

Warnings, Precautions, Contraindications & Adverse Effects

- + Incorrect positioning of the patient on this product can cause:
 - + Pressure marking or ulcers
 - + Neurological damage, which can cause paraplegia
 - + Venous pressure, which can lead to ischaemia
- + Do not immerse in water
- + Do not autoclave
- + Do not use if this product is torn or damaged

Packaging

- + Each Unit is packaged in a recyclable polyethylene bag
- + 1 unit is packaged into a cardboard box for sale

Transport & Storage

- + Transport and store in dry conditions and out of direct light where possible
- + Avoid storing where sharp objects could pierce the skin of the product