




# Technical Data Sheet

## CSM-2550 CS Prone Face Block 5" (Right Intubation Slot)

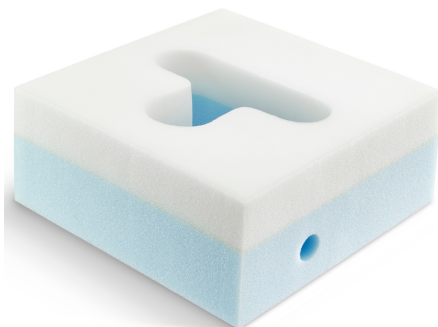

### Description

- Supports patients' faces in the prone position
- Stand alone/ all-in-one product
- Firm base provides stability
- Soft polyurethane foam top for patient comfort and protection
- Distributes the pressure on the patient's cheeks, forehead and chin while in the prone position
- Reduces the incidence of pressure ulcers and marking
- Raises the patient's head to assist with achieving a natural neutral neckline
- Hole in the side of the block to hold tubing
- For patients aged 8 and over, dependent on stature
- For patients weighing up to 200kg

### Application Area

- OR
- Anaesthesia

### Measurements

- Length: 270mm/10.63"
- Width: 246mm/9.69"
- Height: 127mm/5"

### Material Composition

- Top layer of soft polyurethane foam
- Latex & phthalate free

### Order Information

Product Code	Product Name	Qty
CSM-2550	CS Prone Face Block 5" (Right Intubation Slot)	12 Units

### Packaging Information

- Each set is packaged in a recyclable polyethene bag
- 12 units are packaged into a cardboard box for sale

### Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

### Contraindications

Incorrect positioning of the patient on this product can cause:

- Pressure marking or ulcers
- Neurological damage, which can cause paraplegia;
- Venous pressure, which can lead to ischaemia
- Do not use if this product is torn or damaged
- Only to be used by trained medical personnel

### Further Information

- Disposable/Reusable: Disposable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624120084

### Labelling



Manufacturer

UK: Care Surgical Ltd  
Unit 6 Ringtail Road  
Burscough, L40 8JY  
United Kingdom  
Tel: +44 1704 336671

US: Care Surgical LLC  
Marlborough  
MA 01752  
United States  
Tel: +1-866-243-4107

[www.care-surgical.com](http://www.care-surgical.com)



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland

Australian Sponsor, Summit Surgical  
PO Box 60, Landsborough  
QLD 4550 Australia



Medical Device



Do Not Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark

Rev.01 11/09/2023