



Technical Data Sheet

CSM-2715 CS Spine Frame Set



Description

- Supports for patients in the prone position
- Set includes a frame with straps to secure it to the operating table, hip supports and chest support
- Reusable: Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Contains: CS Chest Support, CS Hip Supports (Pair), CS Spine Frame Straps (Set of 4), CS Comfort Covers - 3 Piece (3 sets included)

Application Area

- OR
- Anaesthesia

Measurements

- Length 650mm/25.59"
- Width 510mm/20"
- Height 210mm/8.26"

Material Composition

- Wipeable skin covering
- Polyurethane foam
- Transfer coated polyurethane knitted polyester Acetal plastic
- Carbon Fibre
- Woven Polyester

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM-2715	CS Spine Frame Set	1 Set

Packaging Information

- Each set is packaged in a recyclable polyethene bag
- 1 set is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top
- Avoid storing where sharp objects could pierce the skin of the product

Contraindications

Incorrect positioning of the patient on this product can cause:

- Pressure marking or ulcers
- Neurological damage, which can cause paraplegia
- Venous pressure, which can lead to ischaemia

Examine the pad for damage before and after each use. If the skin is pierced or torn, it cannot be cleaned effectively and should be replaced.

- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel
- If the clips are not fully closed, the support could move during the procedure which could cause injury to the patient or clinician

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624120596

Labelling




Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
Marlborough
MA 01752
United States
Tel: +1-866-243-4107

www.care-surgical.com

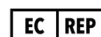

Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland

Australian Sponsor, Summit Surgical
PO Box 60, Landsborough
QLD 4550 Australia



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark

Rev.01 31/08/2023