



# Technical Data Sheet

## CSM-2565 CS Prone Plus Head Support System



### Description

- Reusable
- Helmet holds the CS Prone Plus Face Cushion (CSM-2512) or CS Small Face Cushion (CSM-2514) to maintain their integrity throughout prone procedures
- Space for safe positioning of ET tubes.
- The helmet can be used with or without the mirror
- The mirror enables the clinician to monitor the patient's eyes and tubing during the procedure
- The cogs fitted to the mirror are used to adjust the height or angle of the helmet to achieve the safest natural neutral neckline

### Application Area

- OR
- Anaesthesia
- ICU

### Measurements

- Helmet: L310 | W220 | H90, Mirror: L220 | W220 mm
- Helmet: L12.2 | W8.6 | H3.5, Mirror: L8.6 | W8.6"

### Material Composition

- Helmet made from ABS Plastic
- Acrylic mirror, which incorporates 4 adjustable cogs

### Cleaning

- Clean with standard hospital wipes

### Order Information

Product Code	Product Name	Qty
CSM-2565	CS Prone Plus Head Support System	1 System

### Packaging Information

- Each set is packaged in a recyclable polyethylene bag
- Each set is packaged into a cardboard box for sale

### Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

### Contraindications

Incorrect positioning of the patient on this product can cause:

- POVL
  - pressure ulcers
  - neurological damage
  - venous pressure, which can lead to ischaemia
- 
- Do not use if there is any damage on the helmet or mirror
  - Do not autoclave
  - Only to be used by trained medical personnel

### Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624120121

### Labelling




Manufacturer

UK: Care Surgical Ltd  
Unit 6 Ringtail Road  
Burscough, L40 8JY  
United Kingdom  
Tel: +44 1704 336671

US: Care Surgical LLC  
Marlborough  
MA 01752  
United States  
Tel: +1-866-243-4107

[www.care-surgical.com](http://www.care-surgical.com)

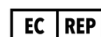

Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland

Australian Sponsor, Summit Surgical  
PO Box 60, Landsborough  
QLD 4550 Australia



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark

Rev.01 7/23