



Technical Data Sheet

CSM-2646 CS Universal Arm Support



Description

- Reusable
- Attaches to most standard and specialist operating tables to support patients' arms in various positions during surgical procedures
- One-step set-up
- 3 rotating joints support full range of movement
- Central clamping mechanism for easy, fast and secure positioning of the arm support
- Including soft foam pad with padded strap
- Spherical joints for smooth rotation

Application Area

- OR
- Anaesthesia

Measurements

- L400 | W150 | H550 mm
- L15.7 | W5.9 | H21.7"

Material Composition

- Stainless steel
- PU-Foam (covered)

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM-2646	CS Universal Arm Support	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag
- Each unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

Incorrect positioning of the patient on this product can cause:

- pressure ulcers
- neurological damage
- venous pressure, which can lead to ischaemia
- Do not use if there is any damage on the product
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624122125

Labelling




Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
Marlborough
MA 01752
United States
Tel: +1-866-243-4107

www.care-surgical.com

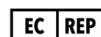

Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland

Australian Sponsor, Summit Surgical
PO Box 60, Landsborough
QLD 4550 Australia



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark

Rev.01 7/23