



Technical Data Sheet

CSM-2699 CS Easy Prone - 150mm - Chest & Pelvis Support



Order Information

Product Code	Product Name	Qty
CSM-2699	CS Easy Prone - 150mm - Chest & Pelvis Support	1 Set

Description

- Set consists of: CSM-2920 CS Easy Prone - 150mm - Chest Only & CSM-2921 CS Easy Prone - 150mm - Pelvis Only
- Raises the patient's chest and pelvis in the prone position
- Enables compression-free positioning of the abdomen
- Reusable, fully radiolucent and fully sealed
- The design prevents pressure in the armpit, neck and genital area
- Suitable for many different body shapes, as it can be flexibly adjusted
- Enables a natural lordosis of the spine
- Allows movements of the table without compromises

Application Area

- OR
- Anaesthesia

Measurements

Chest:

- Length 250mm/9.84"
- Width 500mm/19.69"
- Height 150mm/5.91"

Pelvis:

- Length 336mm/13.23"
- Width 500mm/19.69"
- Height 150mm/5.91"

Material Composition

- Polyurethane foam
- Latex & phthalate free

Packaging Information

- Each component of the set is packaged in a recyclable polyethylene bag
- 1 set is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place

Contraindications

- Incorrect positioning of the patient on the CS Easy Prone can cause:
- Pressure marking or ulcers
- Neurological damage, which can cause paraplegia;
- Venous pressure, which can lead to ischaemia
- Do not use if this product is torn or damaged
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624126994

Cleaning

- Clean using alcohol-based hospital wipes

Labelling




Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
Marlborough
MA 01752
United States
Tel: +1-866-243-4107

www.care-surgical.com

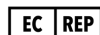

Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland

Australian Sponsor, Summit Surgical
PO Box 60, Landsborough
QLD 4550 Australia



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark

Rev.01 23/11/2023