




Technical Data Sheet

CSM-2913 CS FlexiFin™ Max Stirrups (Pair)

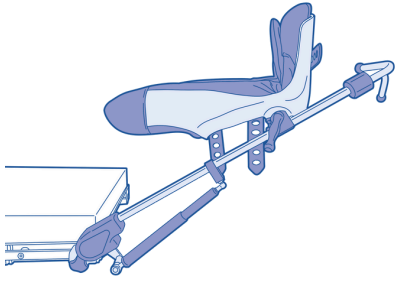

Description

- Heavy stirrup leg holders to protect the patients calves in lithotomy position.
- Includes silicone strap pad attachment system - CS FlexiFin™ Stirrup Pads
- Reusable - clean using alcohol-based hospital wipes
- Latex and phthalate free
- Patient Weight Capacity: 250kg

Application Area

- OR

Measurements

- Length 1000mm/39.37"
- Width 320mm/12.6"
- Height 400mm/15.75"

Material Composition

Stirrup:

- Stainless Steel
- Aluminium
- ABS (Boot)
- Rubber (Shrouds and housings)
- Hook velcro

Pads:

- Polyurethane coated fabric
- Viscoelastic memory foam
- Silicone
- Loop velcro

Cleaning

- Clean using non-alcohol hospital wipes

Order Information

Product Code	Product Name	Qty
CSM-2913	CS FlexiFin™ Max Stirrups (Pair)	1 Pair

Packaging Information

- Each pair is packaged in a recyclable polyethene bag and bubblewrapped
- 1 pair is packaged into a cardboard box for sale

Storage Instructions

- Transport and store in dry conditions

Contraindications

- Incorrect attachment to the surgical tables could result injuries to the patient or nurse
- Do not use if the blade clamp appears to be defective or damaged
- Do not use if the Stirrup Legholders appear to be defective or damaged
- Do not autoclave
- Only to be used by trained medical personnel
- Do not exceed the device safe working load
- Do not immerse in water
- After each use, clean the device with alcohol-based wipes
- Make sure that the device is dry before you store it or use it again

Features

- Patient Weight Capacity: 250kg
- Lithotomy adjustment (from sterile field): Yes
- Lithotomy range (from sterile field): +90° to -45°
- Abduction range (from sterile field): +25° to -10°
- Maximum medial/lateral shift/ boot rotational angle: +-23°
- Product weight: 14,6kg per pair (including pads)

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624129131

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
Marlborough
MA 01752
United States
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Patient Weight Capacity

Rev.01 18/07/2024