



CARE SURGICAL
Specialists in Patient Positioning

Technical Data Sheet

CS Universal Spine Extension™ Set

This Booklet Includes:

EU:

CSM2958 - CS Universal Spine Extension™ SET (EU) Contains:

CSM2955 CS Universal Spine Extension™ - Table Adaptor (EU)

CSM2957 CS Universal Spine Extension™ - Frame

CSM2956 CS Universal Spine Extension™ - Powered Leg

US:

CSM3031 - CS Universal Spine Extension™ SET (US) Contains:

CSM3030 CS Universal Spine Extension™ - Table Adaptor (US)

CSM2957 CS Universal Spine Extension™ - Frame

CSM2956 CS Universal Spine Extension™ - Powered Leg

Asia:

CSM3033 - CS Universal Spine Extension™ SET (Asia) Contains:

CSM3032 CS Universal Spine Extension™ - Table Adaptor (Asia)

CSM2957 CS Universal Spine Extension™ - Frame

CSM2956 CS Universal Spine Extension™ - Powered Leg

AUS:

CSM3035 - CS Universal Spine Extension™ SET (AUS) Contains:

CSM3034 CS Universal Spine Extension™ - Table Adaptor (AUS)

CSM2957 CS Universal Spine Extension™ - Frame

CSM2956 CS Universal Spine Extension™ - Powered Leg

UK:

CSM3037 - CS Universal Spine Extension™ SET (UK) Contains:

CSM3036 CS Universal Spine Extension™ - Table Adaptor (UK)

CSM2957 CS Universal Spine Extension™ - Frame

CSM2956 CS Universal Spine Extension™ - Powered Leg

Japan:

CSM3104 - CS Universal Spine Extension™ SET (Japan) Contains:

CSM3039 CS Universal Spine Extension™ - Table Adaptor (Japan)

CSM2957 CS Universal Spine Extension™ - Frame

CSM2956 CS Universal Spine Extension™ - Powered Leg




CARE SURGICAL
Specialists in Patient Positioning

Technical Data Sheet

CSM2958

CS Universal Spine Extension™ Set

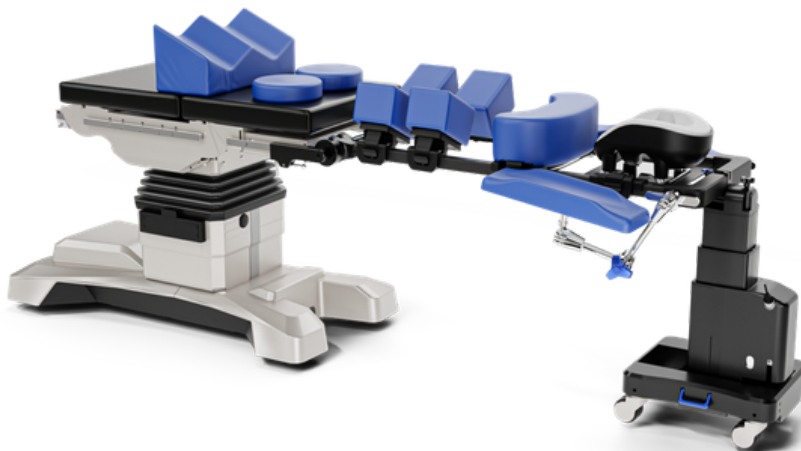



Technical Data Sheet

CSM2958 CS Universal Spine Extension™ Set



Description

- Dual rail extension frame intended to be attached to an operating table, to which specialist supports can be attached to position patients for surgery (e.g. spinal surgery).
- Compatible with sectional pressure management pads and other support attachments (hip/thigh supports).
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1595mm/62.78"
- Width: 406mm/15.98"
- Height: 960mm/37.80"(Max) | 560mm/22.05" (Min)

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium
- Polypropylene
- Rubber
- Impact resistant plastic

Power Supply

- Power Supply - 100-240 VAC, 50/60 Hz
- Standby Power - Approx. 0.5 W (120/230 V)

Order Information

Product Code	Product Name	Qty
CSM2958	CS Universal Spine Extension™ Set	1 Set

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the extension frame is securely attached to the table before use
- Check if the extension frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624126000

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM2956 CS Universal Spine Extension™ - Powered Leg



Description

- Electrical telescopic leg
- Moves respective to the surgical table
- used with CSM3030 and CSM2957
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 460mm/18.11"
- Width: 397mm/15.59"
- Height: 1030mm/40.55"(Max) | 630mm/24.80" (Min)

Material Composition

- Stainless steel
- Polypropylene
- Rubber
- Aluminium
- Impact resistant plastic

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2956	CS Universal Spine Extension™ - Powered Leg	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure brakes are engaged before product is in use
- Check if the leg is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624126017

Power Supply

- Power Supply - 100-240 VAC, 50/60 Hz
- Standby Power - Approx. 0.5 W (120/230 V)

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM2957 CS Universal Spine Extension™ - Frame



Description

- Frame supports the patient and supports positioning equipment
- Connects to CSM3030 and CSM2956
- Clamp positioning equipment onto frame members
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1483mm/58.39"
- Width: 503mm/19.80"
- Height: 78mm/3.07"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2957	CS Universal Spine Extension™ - Frame	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the frame is securely attached to the table adaptor before use
- Check if the frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624125997

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Technical Data Sheet

CSM2955 CS Universal Spine Extension™ - Table Adaptor EU



Order Information

Product Code	Product Name	Qty
CSM2955	CS Universal Spine Extension™ - Table Adaptor EU	1 Unit

Description

- Adaptable rail clamp that supports different width surgical tables
- Clamp onto surgical table rails. To be used with CSM2956 and CSM2957. CSM2957 (frame) attaches to CSM2955 (adaptor).
- For use on compatible surgical tables in surgical environments and ICUs
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 226mm/8.1"
- Width: 740.5mm/29.07"(Max) | 540.5mm/21.2" (Min)
- Height: 170mm/7.10"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the adaptor is securely attached to the table before use
- Check the adaptor for damage before and after each use. Do not use if damaged Contact supplier for advice
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- UDI-DI/GTIN: 5060624129551

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark




CARE SURGICAL
Specialists in Patient Positioning

Technical Data Sheet

CSM3031

CS Universal Spine Extension™ Set

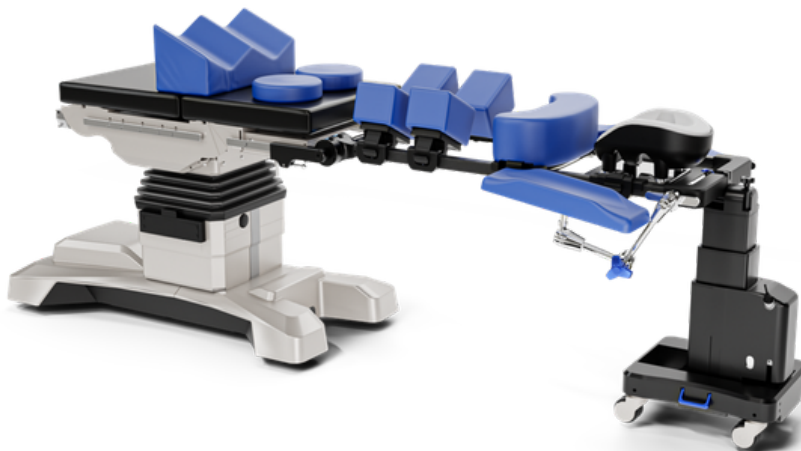



Technical Data Sheet

CSM3031 CS Universal Spine Extension™ Set (US)



Description

- Dual rail extension frame intended to be attached to an operating table, to which specialist supports can be attached to position patients for surgery (e.g. spinal surgery).
- Compatible with sectional pressure management pads and other support attachments (hip/thigh supports).
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1595mm/62.78"
- Width: 406mm/15.98"
- Height: 960mm/37.80"(Max) | 560mm/22.05" (Min)

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium
- Polypropylene
- Rubber
- Impact resistant plastic

Power Supply

- Power Supply - 100-240 VAC, 50/60 Hz
- Standby Power - Approx. 0.5 W (120/230 V)

Order Information

Product Code	Product Name	Qty
CSM3031	CS Universal Spine Extension™ Set (US)	1 Set

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the extension frame is securely attached to the table before use
- Check if the extension frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624126536

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Technical Data Sheet

CSM2956 CS Universal Spine Extension™ - Powered Leg



Description

- Electrical telescopic leg
- Moves respective to the surgical table
- used with CSM3030 and CSM2957
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 460mm/18.11"
- Width: 397mm/15.59"
- Height: 1030mm/40.55"(Max) | 630mm/24.80" (Min)

Material Composition

- Stainless steel
- Polypropylene
- Rubber
- Aluminium
- Impact resistant plastic

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2956	CS Universal Spine Extension™ - Powered Leg	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure brakes are engaged before product is in use
- Check if the leg is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624126017

Power Supply

- Power Supply - 100-240 VAC, 50/60 Hz
- Standby Power - Approx. 0.5 W (120/230 V)

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Technical Data Sheet

CSM2957 CS Universal Spine Extension™ - Frame



Description

- Frame supports the patient and supports positioning equipment
- Connects to CSM3030 and CSM2956
- Clamp positioning equipment onto frame members
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1483mm/58.39"
- Width: 503mm/19.80"
- Height: 78mm/3.07"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2957	CS Universal Spine Extension™ - Frame	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the frame is securely attached to the table adaptor before use
- Check if the frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624125997

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Technical Data Sheet

CSM3030 CS Universal Spine Extension™ - Table Adaptor US



Order Information

Product Code	Product Name	Qty
CSM3030	CS Universal Spine Extension™ - Table Adaptor US	1 Unit

Description

- Adaptable rail clamp that supports different width surgical tables
- Clamp onto surgical table rails. To be used with CSM2956 and CSM2957. CSM2957 (frame) attaches to CSM3030 (adaptor).
- For use on compatible surgical tables in surgical environments and ICUs
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 226mm/8.1"
- Width: 740.5mm/29.07"(Max) | 540.5mm/21.2" (Min)
- Height: 170mm/7.10"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the adaptor is securely attached to the table before use
- Check the adaptor for damage before and after each use. Do not use if damaged Contact supplier for advice
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624129551

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark




CARE SURGICAL
Specialists in Patient Positioning

Technical Data Sheet

CSM3033

CS Universal Spine Extension™ Set

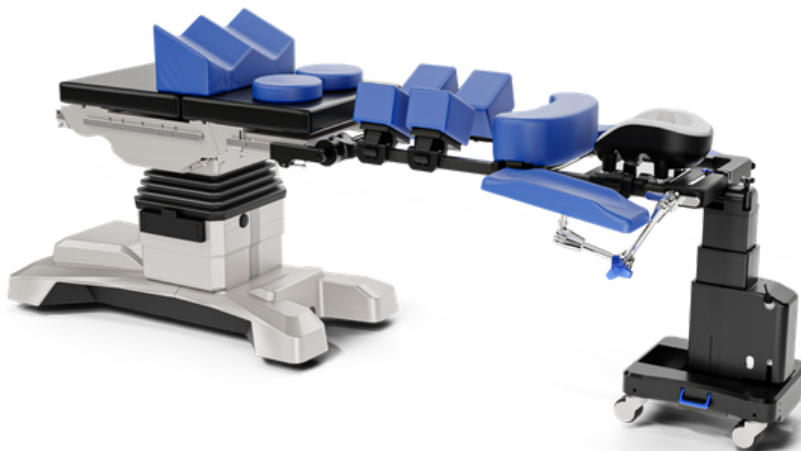



Technical Data Sheet

CSM3033 CS Universal Spine Extension™ Set (Asia)



Description

- Dual rail extension frame intended to be attached to an operating table, to which specialist supports can be attached to position patients for surgery (e.g. spinal surgery).
- Compatible with sectional pressure management pads and other support attachments (hip/thigh supports).
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Max Dims - L 1595mm | W 406mm | H 960mm
- Min Dims- L1595mm | W 406mm | H 560mm

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium
- Polypropylene
- Rubber
- Impact resistant plastic

Cleaning

- Clean with standard hospital wipes

Power Supply

- Power Supply - 100-240 VAC, 50/60 Hz
- Standby Power - Approx. 0.5 W (120/230 V)

Order Information

Product Code	Product Name	Qty
CSM3033	CS Universal Spine Extension™ Set (Asia)	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the extension frame is securely attached to the table before use
- Check if the extension frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624126550

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM2956 CS Universal Spine Extension™ - Powered Leg



Description

- Electrical telescopic leg
- Moves respective to the surgical table
- used with CSM 2955 and CSM 2957
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Max Dims: L 960mm | W 460mm | H 406mm
- Min Dims: L 560mm | W 460mm | H 406mm

Material Composition

- Stainless steel
- Polypropylene
- Rubber
- Aluminium
- Impact resistant plastic

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2956	CS Universal Spine Extension™ - Powered Leg	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure brakes are engaged before product is in use
- Check if the leg is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624126017

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM2957 CS Universal Spine Extension™ - Frame



Description

- Frame supports the patient and supports positioning equipment
- Connects to CSM 2955 and CSM 2956
- Clamp positioning equipment onto frame members
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

L 1482.5mm | W 502.4mm | H 77.61mm

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2957	CS Universal Spine Extension™ - Frame	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the frame is securely attached to the table adaptor before use
- Check if the frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624125997

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Technical Data Sheet

CSM3032 CS Universal Spine Extension™ - Table Adaptor (Asia)



Order Information

Product Code	Product Name	Qty
CSM3032	CS Universal Spine Extension™ - Table Adaptor (Asia)	1 Unit

Description

- Adaptable rail clamp that supports different width surgical tables
- Clamp onto surgical table rails. To be used with CSM2956 and CSM2957. CSM2957 (frame) attaches to CSM3032 (adaptor).
- For use on compatible surgical tables in surgical environments and ICUs
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Max Dims: L 226mm | W 540.5/740.5mm | H 150.5mm

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the adaptor is securely attached to the table before use
- Check the adaptor for damage before and after each use. Do not use if damaged Contact supplier for advice
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- Basic UDI/GTIN No: 5060624126543

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark




CARE SURGICAL
Specialists in Patient Positioning

Technical Data Sheet

CSM3035

CS Universal Spine Extension™ Set

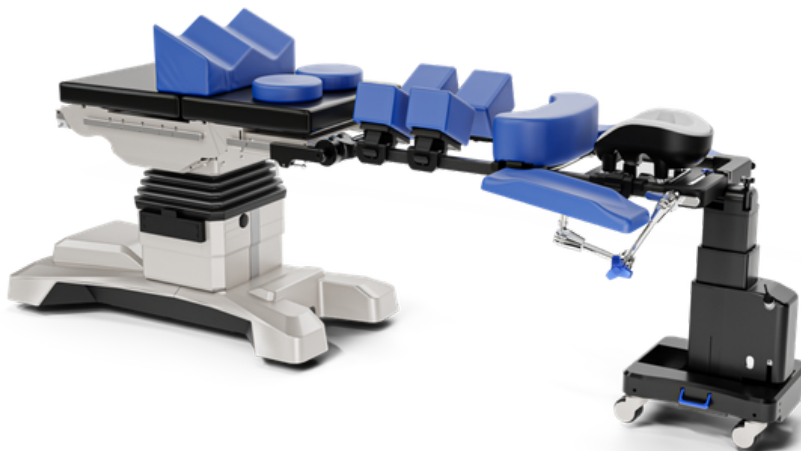



Technical Data Sheet

CSM3035 CS Universal Spine Extension™ Set (AUS)



Description

- Dual rail extension frame intended to be attached to an operating table, to which specialist supports can be attached to position patients for surgery (e.g. spinal surgery).
- Compatible with sectional pressure management pads and other support attachments (hip/thigh supports).
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1595mm/62.78"
- Width: 406mm/15.98"
- Height: 960mm/37.80"(Max) | 560mm/22.05" (Min)

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium
- Polypropylene
- Rubber
- Impact resistant plastic

Power Supply

- Power Supply - 100-240 VAC, 50/60 Hz
- Standby Power - Approx. 0.5 W (120/230 V)

Order Information

Product Code	Product Name	Qty
CSM3035	CS Universal Spine Extension™ Set (AUS)	1 Unit

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the extension frame is securely attached to the table before use
- Check if the extension frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN 5060624127328

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland

Australian Sponsor:

MDSS AU Pty Ltd, Ground Floor, Suite 55,
97-99 Bathurst Street, Sydney, NSW 2000, Australia



Technical Data Sheet

CSM2956 CS Universal Spine Extension™ - Powered Leg



Description

- Electrical telescopic leg
- Moves respective to the surgical table
- used with CSM3030 and CSM2957
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 460mm/18.11"
- Width: 397mm/15.59"
- Height: 1030mm/40.55"(Max) | 630mm/24.80" (Min)

Material Composition

- Stainless steel
- Polypropylene
- Rubber
- Aluminium
- Impact resistant plastic

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2956	CS Universal Spine Extension™ - Powered Leg	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure brakes are engaged before product is in use
- Check if the leg is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- UDI-DI/GTIN 5060624126017

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Technical Data Sheet

CSM2957 CS Universal Spine Extension™ - Frame



Description

- Frame supports the patient and supports positioning equipment
- Connects to CSM3030 and CSM2956
- Clamp positioning equipment onto frame members
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1483mm/58.39"
- Width: 503mm/19.80"
- Height: 78mm/3.07"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2957	CS Universal Spine Extension™ - Frame	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the frame is securely attached to the table adaptor before use
- Check if the frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN 5060624125997

Power Supply

- Power Supply - 100-240 VAC, 50/60 Hz
- Standby Power - Approx. 0.5 W (120/230 V)

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Technical Data Sheet

CSM3034 CS Universal Spine Extension™ - Table Adaptor AUS



Order Information

Product Code	Product Name	Qty
CSM3034	CS Universal Spine Extension™ - Table Adaptor AUS	1 Unit

Description

- Adaptable rail clamp that supports different width surgical tables
- Clamp onto surgical table rails. To be used with CSM2956 and CSM2957. CSM2957 (frame) attaches to CSM3034 (adaptor).
- For use on compatible surgical tables in surgical environments and ICUs
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 226mm/8.1"
- Width: 740.5mm/29.07"(Max) | 540.5mm/21.2" (Min)
- Height: 170mm/7.10"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the adaptor is securely attached to the table before use
- Check the adaptor for damage before and after each use. Do not use if damaged Contact supplier for advice
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- UDI-DI/GTIN: 5060624128158

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland

Australian Sponsor:
MDSS AU Pty Ltd, Ground Floor, Suite 55,
97-99 Bathurst Street, Sydney, NSW 2000, Australia



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark




CARE SURGICAL
Specialists in Patient Positioning

Technical Data Sheet

CSM3037

CS Universal Spine Extension™ Set

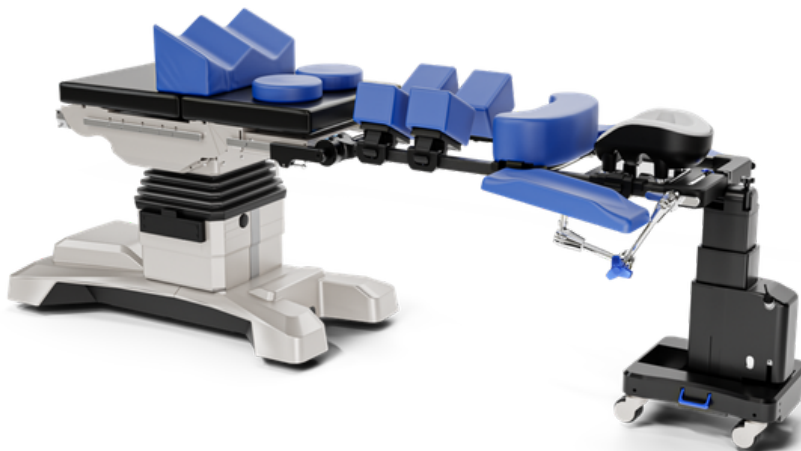



Technical Data Sheet

CSM3037 CS Universal Spine Extension™ Set (UK)



Description

- Dual rail extension frame intended to be attached to an operating table, to which specialist supports can be attached to position patients for surgery (e.g. spinal surgery).
- Compatible with sectional pressure management pads and other support attachments (hip/thigh supports).
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1595mm/62.78"
- Width: 406mm/15.98"
- Height: 960mm/37.80"(Max) | 560mm/22.05" (Min)

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium
- Polypropylene
- Rubber
- Impact resistant plastic

Order Information

Product Code	Product Name	Qty
CSM3037	CS Universal Spine Extension™ Set (UK)	1 Unit

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the extension frame is securely attached to the table before use
- Check if the extension frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624127335

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark



Technical Data Sheet

CSM2956 CS Universal Spine Extension™ - Powered Leg



Description

- Electrical telescopic leg
- Moves respective to the surgical table
- used with CSM3030 and CSM2957
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 460mm/18.11"
- Width: 397mm/15.59"
- Height: 1030mm/40.55"(Max) | 630mm/24.80" (Min)

Material Composition

- Stainless steel
- Polypropylene
- Rubber
- Aluminium
- Impact resistant plastic

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2956	CS Universal Spine Extension™ - Powered Leg	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure brakes are engaged before product is in use
- Check if the leg is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624126017

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM2957 CS Universal Spine Extension™ - Frame



Description

- Frame supports the patient and supports positioning equipment
- Connects to CSM3030 and CSM2956
- Clamp positioning equipment onto frame members
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1483mm/58.39"
- Width: 503mm/19.80"
- Height: 78mm/3.07"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2957	CS Universal Spine Extension™ - Frame	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the frame is securely attached to the table adaptor before use
 - Check if the frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
 - Do not exceed the maximum weight mentioned on the IFU
 - Do not immerse in water
 - Do not autoclave
 - Only to be used by trained medical personnel
- L 1482.5mm | W 502.4mm | H 77.61mm

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- UDI-DI/GTIN: 5060624125997

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM3036 CS Universal Spine Extension™ - Table Adaptor UK



Description

- Adaptable rail clamp that supports different width surgical tables
- Clamp onto surgical table rails. To be used with CSM2956 and CSM2957. CSM2957 (frame) attaches to CSM3036 (adaptor).
- For use on compatible surgical tables in surgical environments and ICUs
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 226mm/8.1"
- Width: 738.5mm/29.07"(Max) | 538.5mm/21.28" (Min)
- Height: 180.5mm/7.10"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Order Information

Product Code	Product Name	Qty
CSM3036	CS Universal Spine Extension™ - Table Adaptor UK	1 Unit

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the adaptor is securely attached to the table before use
- Check the adaptor for damage before and after each use. Do not use if damaged Contact supplier for advice
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- UDI-DI/GTIN: 5060624128165

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark




CARE SURGICAL
Specialists in Patient Positioning

Technical Data Sheet

CSM3104

CS Universal Spine Extension™ Set

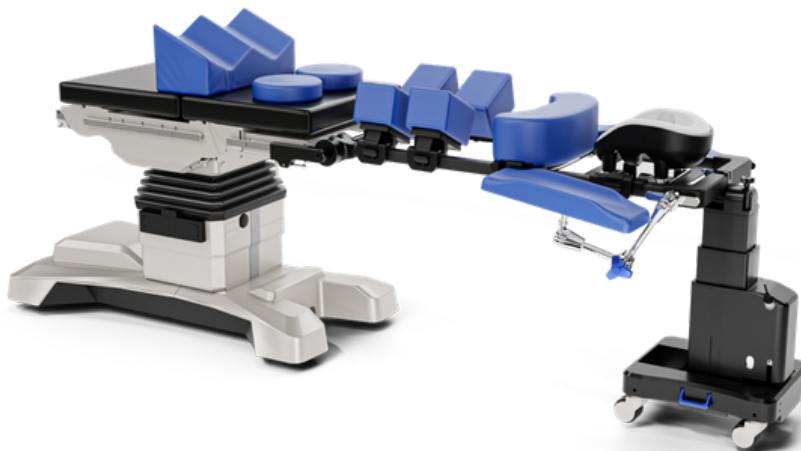



Technical Data Sheet

CSM3104 CS Universal Spine Extension™ SET (Japan)



Description

- Dual rail extension frame intended to be attached to an operating table, to which specialist supports can be attached to position patients for surgery (e.g. spinal surgery).
- Compatible with sectional pressure management pads and other support attachments (hip/thigh supports).
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1595mm/62.78"
- Width: 406mm/15.98"
- Height: 960mm/37.80"(Max) | 560mm/22.05" (Min)

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium
- Polypropylene
- Rubber
- Impact resistant plastic

Order Information

Product Code	Product Name	Qty
CSM3104	CS Universal Spine Extension™ Set (Japan)	1 Unit

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the extension frame is securely attached to the table before use
- Check if the extension frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624127335

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM2956 CS Universal Spine Extension™ - Powered Leg



Description

- Electrical telescopic leg
- Moves respective to the surgical table
- used with CSM3030 and CSM2957
- Clean using non-alcohol hospital wipes
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 460mm/18.11"
- Width: 397mm/15.59"
- Height: 1030mm/40.55"(Max) | 630mm/24.80" (Min)

Material Composition

- Stainless steel
- Polypropylene
- Rubber
- Aluminium
- Impact resistant plastic

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2956	CS Universal Spine Extension™ - Powered Leg	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure brakes are engaged before product is in use
- Check if the leg is damaged before and after each use. Do not use if damaged and contact supplier for advice
- Do not exceed the maximum weight mentioned on the IFU
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA: Y
- MD Class: 1
- UDI-DI/GTIN: 5060624126017

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM2957 CS Universal Spine Extension™ - Frame



Description

- Frame supports the patient and supports positioning equipment
- Connects to CSM3030 and CSM2956
- Clamp positioning equipment onto frame members
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 1483mm/58.39"
- Width: 503mm/19.80"
- Height: 78mm/3.07"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Cleaning

- Clean with standard hospital wipes

Order Information

Product Code	Product Name	Qty
CSM2957	CS Universal Spine Extension™ - Frame	1 Unit

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the frame is securely attached to the table adaptor before use
 - Check if the frame is damaged before and after each use. Do not use if damaged and contact supplier for advice
 - Do not exceed the maximum weight mentioned on the IFU
 - Do not immerse in water
 - Do not autoclave
 - Only to be used by trained medical personnel
- L 1482.5mm | W 502.4mm | H 77.61mm

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- UDI-DI/GTIN: 5060624125997

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity
Assessed



CE Mark



Technical Data Sheet

CSM3039 CS Universal Spine Extension™ - Table Adaptor (Japan)



Description

- Adaptable rail clamp that supports different width surgical tables
- Clamp onto surgical table rails. To be used with CSM2956 and CSM2957. CSM2957 (frame) attaches to CSM3036 (adaptor).
- For use on compatible surgical tables in surgical environments and ICUs
- Clean using non-alcohol hospital wipes
- Moves respective to the surgical table
- Number of uses: Unrestricted, subject to damage (see damage indicators in the Contraindications column)
- Latex and phthalate free
- Instructions for use are packaged with the product

Application Area

- OR
- Anaesthesia

Measurements

- Length: 226mm/8.1"
- Width: 738.5mm/29.07"(Max) | 538.5mm/21.28" (Min)
- Height: 180.5mm/7.10"

Material Composition

- Stainless steel
- Carbon Fibre
- Aluminium

Order Information

Product Code	Product Name	Qty
CSM3039	CS Universal Spine Extension™ - Table Adaptor (Japan)	1 Unit

Cleaning

- Clean with standard hospital wipes

Packaging Information

- Each unit is packaged in a recyclable polyethylene bag and bubblewrap
- 1 unit is packaged into a cardboard box for sale

Storage Instructions

- First in first out
- Store in a cool dry place
- Avoid direct sunlight
- Do not store anything heavy on top

Contraindications

- Ensure the adaptor is securely attached to the table before use
- Check the adaptor for damage before and after each use. Do not use if damaged Contact supplier for advice
- Do not immerse in water
- Do not autoclave
- Only to be used by trained medical personnel

Further Information

- Disposable/Reusable: Reusable
- Sterility: N
- CE/FDA/TGA/UKCA:Y
- MD Class: 1
- UDI-DI/GTIN: 5060624128165

Labelling



Manufacturer

UK: Care Surgical Ltd
Unit 6 Ringtail Road
Burscough, L40 8JY
United Kingdom
Tel: +44 1704 336671

US: Care Surgical LLC
43 Broad Street
Suite B106, Hudson
MA 01749, USA
Tel: +1-866-243-4107

www.care-surgical.com



Care Surgical Supplies Ltd, Office 2, 12a Lower Main Street, Lucan Co. Dublin, K78 X5P8 Ireland



Medical Device



Reuse



Latex Free



Non Sterile



UK Conformity Assessed



CE Mark